USER MANUAL

Mat - massage mattress Malatec 23723

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Mat – Massagematratze Malatec 23723

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Mat - masážní matrace Malatec 23723

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Mat - matelas de massage Malatec 23723

INSTRUCCIONES DE OPERACION

Mat - colchón de masaje Malatec 23723

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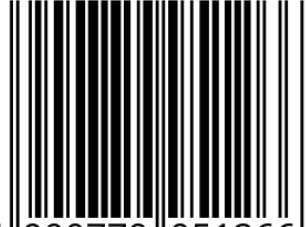
5 900779 951266 >

EN

Ladies and Gentlemen, thank you for purchasing our product!
Before starting to use, please read the following instructions
for the correct use of the product. Please keep these instructions on
future and to follow its recommendations, as non-compliance may endanger life or health.

Important: KEEP FOR FUTURE REFERENCE!

- Please read it completely before use and follow all instructions in the manual.
- Do not use the product if any part is missing or damaged.
- Please keep this manual for future reference.

SPECIFICATION:

170 x 60 x 3 cm
Input: 100-240V AC, 50/60 Hz

CHARACTERISTIC

- The set includes: a massage mattress with a controller (remote control), an adapter
- With a smart, elegant design, the massage mattress is equipped with high-frequency massage balls to help relieve body aches, reduce stress, discomfort, fatigue, help you relax your body and mind quickly.
- Intelligent protection mode helps the mattress automatically disconnect the power supply after 15 minutes of continuous use to prevent the motor from overheating.
- Durable, strong material is comfortable to use.
- The mattress is compact, lightweight, easy to carry and can be used on chairs, beds, sofas.
- The product is suitable for people who often suffer from low back pain, the elderly and those who work a lot.
- Regular massage supports the activity of immune cells and brain activity. Thus, it supports blood circulation, increases oxygenation of cells and improves sleep.

USE:

- 1 Place the massage cushion on a chair, bed, sofa or floor.
- 2 Connect the adapter cord to the remote control.
- 3 Plug the adapter into an electrical outlet, press the main switch of the remote control to start and select massage modes and intensity.
- 4 When you stop using it, turn off the switch on the remote control and disconnect the power.
- 5 Allow the mattress to cool down before storing or cleaning.
- 6 The mattress should be cleaned with a dry cloth only, and before cleaning, disconnect the power from the socket.
- 7 Store the mattress in a cool, dry place.

ATTENTION

- The product has an automatic shutdown mode after 15 minutes of operation to prevent overheating. The product must cool down before using it again. Continuous use may cause the product to overheat, resulting in damage or shortening product life.
- Must be supervised when used by children, sick people and people with disabilities
- People using a pacemaker should consult a physician before using the product.
- Use only the original adapter that came with the product.
- Do not connect the product when not in use. If the product is not in use, disconnect it.
- Always disconnect the power cord after use or before cleaning.
- Do not use or store the product in a place where it might get wet.
- Do not grab the power cord and pull the product.
- Do not place the power cord near a heat source.
- Do not use outdoors.
- Do not drop or put anything in the mattress.
- Do not use the product near places where you use flammable aerosols and sprays.
- Do not run and jump on the product.
- Cases where serious injury may occur occur when instructions are ignored.

WARNINGS

- 1 Do not use the product directly on exposed skin to prevent burns, abrasions, etc.
- 2 When you feel pain or discomfort, stop using the product immediately
- 3 Do not use the product with injuries to the body and in people with serious diseases, such as malignant tumors, high fever and skin diseases, consult a doctor before using the preparation.
- 4 Minors should use the product under the supervision of parents or adults who may direct the use of the product.

Information for the consumer on the handling of waste electrical and electronic equipment:

Waste electrical and electronic equipment may not be placed with other household waste. Such equipment should be collected selectively, which is indicated by the attached marking - a crossed-out wheeled waste container. Failure to comply with the above rule, in the case of improper disposal of used equipment, may pose a threat to the natural environment and human health, resulting from the presence of hazardous components in the equipment, such as: electrical wiring, plastic, switches, etc. In order to avoid the above threat, these components should be collected and properly processed by specialized companies. Self-disassembly of worn equipment is unacceptable. Households play an important role in contributing to reuse and recovery, including recycling of waste equipment. This is done in particular through active participation in the waste equipment collection system.

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